




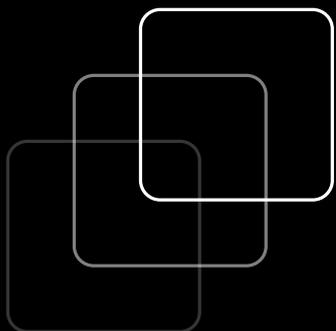

RUBBER SHEET COLLECTION

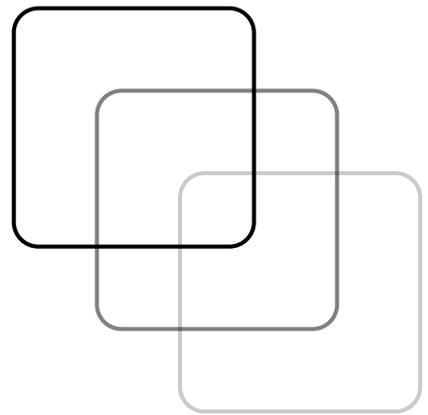

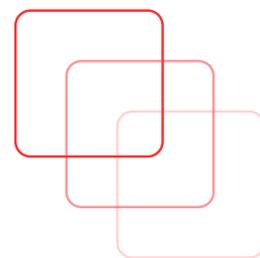

About Us

Shikori Group is a leading manufacturer and supplier of premium mackintosh rubber sheets, widely known as hospital rubber sheets. Our products are designed to provide hygienic, waterproof mattress protection and durable solutions for hospitals, clinics, nursing homes, and home care environments. With years of experience, we combine advanced manufacturing processes with high quality rubber fabric to ensure consistent performance and safety.

Our commitment to excellence extends to timely delivery, customizable options, and dedicated customer support. Every rubber sheet undergoes rigorous quality checks to meet healthcare standards, ensuring reliability, comfort, and long lasting use. **Shikori Group** is proud to be a trusted rubber sheet supplier for medical institutions, distributors, and care facilities across India and international markets.

SHIKORI
GROUP


Product Portfolio



Mackintosh Rubber Sheet **(Hospital rubber sheet)**

Shikori's Hospital Rubber Sheet (Mackintosh sheet) is bright colored, high quality flexible fabric, rubber laminated from both sides, available in one side Blue/Green, other side Red in our standard range.

Manufacturing Plants Photos

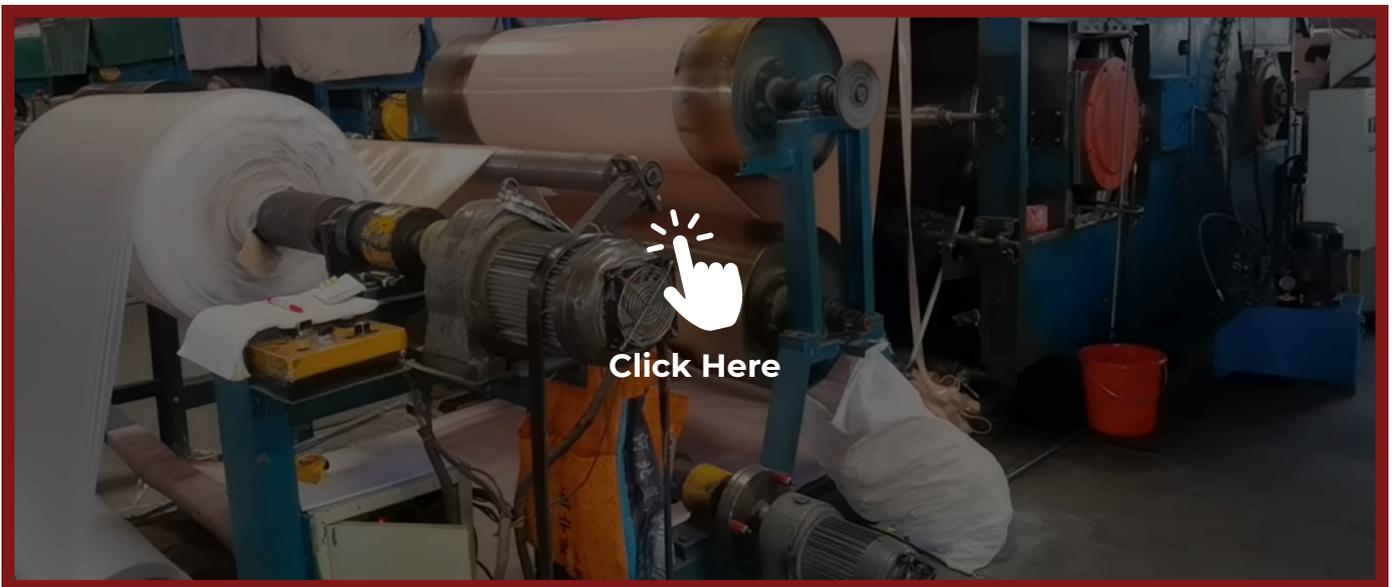

Key Features



100% Waterproof & Leak-proof



High Tensile Strength & Tear-Resistant



Made from Premium Rubber & Fabric Lamination



Easy to Clean | Reusable | Washable



Anti-Microbial & Anti-Allergic Options Available



Temperature Resistant & Flexible



Available in Rolls & Pre-Cut Sizes




Long-life durability for repeated use




Suitable for Autoclaving (optional grade)

Benefit

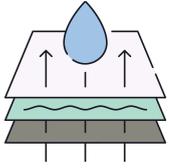

Complete Waterproof Protection

Prevents liquid penetration and protects mattresses from spills, stains, and contamination, ensuring a hygienic patient environment.

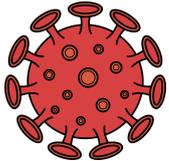

Enhanced Infection Control

Acts as a protective barrier against fluids, helping maintain cleanliness and supporting hospital hygiene standards.



Durable & Long-Lasting

Manufactured with high-quality rubber fabric that withstands regular hospital use, repeated cleaning, and daily wear.



Reusable & Cost-Effective

Easy to clean and disinfect, making it a sustainable and economical solution for hospitals and healthcare facilities.




Flexible & Comfortable

Soft, flexible material adapts to patient beds without discomfort while maintaining protective performance.



Multi-Application Use

Suitable for hospital beds, examination tables, nursing homes, clinics, and home care settings.

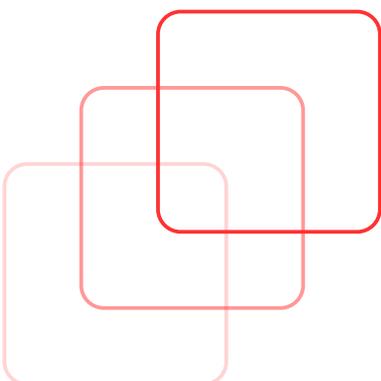

Real Business Impact



Infection Control & Patient Safety

Helps hospitals maintain strict hygiene standards with waterproof and easy-to-sanitize surfaces.



Cost Reduction for Healthcare Facilities

Durable rubber sheets minimize frequent replacements, lowering long-term operational costs.



Increased Operational Efficiency

Easy cleaning and quick drying save staff time and improve turnaround of hospital beds.



Strong Dealer & Distributor Margins

High demand in hospitals, nursing homes, and home care creates consistent repeat orders.



Compliance & Quality Assurance

Manufactured under strict quality checks aligned with healthcare standards.



Scalability for Large Institutions

Capable of fulfilling bulk orders with consistent quality and timely delivery.



Brand Trust & Long-Term Partnerships

Positioning as a dependable supplier ensures recurring institutional contracts.

Our Range



SPECIFICATION

Standard Grade

Range	Width	Colour
Standard-RB	90 CM	Red x Blue
Standard-RG	90 CM	Red x Green

Premium Grade

Range	Width	Colour
Premium-RB	90 CM	Red x Blue
Premium-RG	90 CM	Red x Green
Premium-RB-(Wide)	110 CM	Red x Blue

Deluxe Grade

Range	Width	Colour	Autoclavable
Deluxe-Supreme	90 CM	Red x Blue	Yes
Deluxe-Ultra	90 CM	Red x Green	Yes
Min. 5 Times Autoclavable			

Format: Available in piece and roll form

Roll length: 10–35 meters (depending on variant)

Customisation options are available on request.

Inquiry Now



Applications



Hospital & Clinic Beds



ICU & Maternity Units



Pediatric & Neonatal Units



Examination Table Covers



Nursing Homes




Physiotherapy centers



Home medical care / Senior care

Here is a comprehensive list of all possible applications of Hospital Rubber sheets (Mackintosh sheets) across different industries and usage segments.

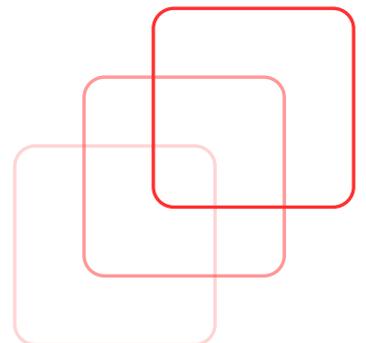

Applications of Hospital Mackintosh Rubber sheets



Healthcare & Medical Sector

- Hospital beds mattress protector
- ICU beds waterproof protection
- Operation theatre table covering
- Labour & maternity room beds
- Delivery beds in gynecology departments
- Neonatal / Pediatric bed protection
- Patient examination tables
- Emergency room & trauma beds
- Physiotherapy & orthopedic tables
- Nursing homes & old-age homes
- Burn unit beds (fluid containment)
- Post-surgery patient care beds
- Cleaning room surface protection
- Stretchers & ambulance beds
- Dental clinic chair covers
- Veterinary hospitals & surgeries

Home & Personal Healthcare

- Home nursing bed protectors
- Senior citizen patient care beds
- Incontinence support beds
- Kids, infant & baby bed protection
- Protective underlay for bedridden patients

Institutional & Commercial Use

- Schools & hostel sick rooms
- Physiotherapy & sports injury centers
- Ayurvedic massage centers
- Spa & massage table covers
- Beauty salon waxing tables
- Tattoo studio seating protection
- Gym physiotherapy and recovery tables

Veterinary & Animal Care

- Veterinary hospitals & clinics
- Animal surgery tables
- Pet care grooming tables
- Cattle farm birthing and treatment beds

Product Care Instructions

Handling

- Unroll the roll and cut the piece into desired length.
- Clean the applied powder by wiping out with soft cotton cloth.
- Spread a hospital rubber sheet (Mackintosh) on the bed with desired colour upside.
- It is recommended to put a cotton bed sheet over a hospital rubber sheet (Mackintosh).
- Keep away any sharp edge material from the hospital rubber sheet (Mackintosh).
- Avoid prolonged contact from any type of oil/waxes.

Disinfection

- In general, Shikori's rubber coating resists all kinds of disinfectants.
- Generally surface wiping can be done numerous times without any problems.
- We advise you to test a particular surface disinfectant in advance or contact us in case of any doubts.

Sterilization/Autoclavation

- Shikori's coatings have the best possible resistance against autoclave sterilization.
- All qualities have been tested and approved for autoclave cycles. Limits can be advised as per quality.

Washing

- Hospital rubber sheets (Mackintosh) can be easily washed for > 100 times as and when required with mild soap/detergent in cold water (excluding possible physical damage of the coating due to contact with sharp abrasive surfaces while washing).
- Washing must be done with neutral detergents; softening agents are not allowed.
- Machine wash or brushing is not permitted.

Drying

Adequate drying is a crucial factor in the care and maintenance of coated articles.

- Wipe dry after washing.
- Drying may be affected by hanging out.
- Do not press, mangle, spin or tumble.
- Do not dry in direct sunlight.
- No ironing.

Storage/preservation

- When not in use, sprinkle or apply talcum powder and keep in roll form duly wrapped in paper, polythene cover.
- To be kept in a cool and dry place, free from ozone with temperature ranging between 10 °C to 35 °C.
- Avoid direct sunlight, UV radiation and contact with hydrocarbon, solvent, acid-alkalis for a prolonged period.

Handling

Surface Preparation

Clean the rubber sheet regularly using mild soap and lukewarm water. Wipe gently with a soft cloth and allow it to dry completely before reuse. Avoid harsh chemicals, bleach, or abrasive cleaners to maintain surface quality and waterproof performance.

Protective Layer

For enhanced comfort and hygiene, it is recommended to place a cotton bed sheet over the rubber sheet during use.

Avoid Sharp Objects

Keep the sheet away from sharp edges, pointed instruments, or rough surfaces that may cause punctures or tears.



Storage

Clean Before Storage

Ensure the rubber sheet is completely clean and dry before storing. Moisture trapped during storage may affect hygiene and material performance.

Keep Away from Chemicals


Store away from oils, solvents, disinfectants, and corrosive substances that may damage the rubber surface.

Avoid Tight Folding

Do not fold tightly for long periods. Roll the sheet loosely to prevent surface creases, cracks, or stress marks.



Installation Procedure

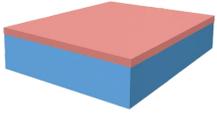

Surface Preparation

Ensure the mattress or examination table is clean, dry, and free from dust or sharp particles before installation.



Unroll & Position

Unroll the mackintosh rubber sheet carefully and place it evenly over the mattress with the desired color side facing upward.



Proper Alignment

Adjust the sheet so it covers the required area fully. Make sure there are no wrinkles or folds to ensure smooth and comfortable placement.



Secure the Edges

Tuck the edges firmly under the mattress or bed frame to prevent movement during patient use.



Add Top Layer (Recommended)

For enhanced comfort and hygiene, place a cotton bedsheet over the rubber sheet before patient use.



Final Check

Inspect the installation to confirm proper coverage, smooth alignment, and secure positioning before use.

Pallet Packing & Container Loading

At Shikori, we ensure that every product is packed with precision and care to maintain its quality during transit. Our export packing is designed to withstand long-distance transportation and harsh shipping conditions.

Pallet Packing Features:

- Durable Material
- Protective Wrapping
- Weight Optimization

Container Loading Options:

- Strategic Arrangement
- Quality Checks
- Secure Loading Techniques



Certificate



Frequently Asked Questions

1

What is a Mackintosh Hospital Rubber Sheet?

A waterproof laminated rubber sheet used to protect hospital mattresses and beds from fluid contamination while maintaining hygiene and durability.

2

Do Shikori Mackintosh sheets meet hospital quality and international standards

Yes. Products comply with ISO 9001, CE requirements, and quality testing is conducted for thickness, elasticity, strength, and waterproofing. Third-party test reports available on request.

3

What sizes and thickness options are available?

Available in 90–150 cm widths, 10–30 meter rolls, and 0.40–1.0 mm thickness with custom sizing and OEM options.

4

Is the sheet washable and reusable?

Yes, it is fully washable, reusable, tear-resistant, and designed for repeated hospital-grade cleaning without losing waterproof performance.

5

What is the minimum order quantity?

MOQ starts from 500 meters, flexible for distributors, export buyers, and new market entry requirements.



SHIKORI
GROUP

SHIKORI GROUP OF COMPANIES



+91 9403268796



export@shikorigroup.com



www.shikorigroup.com



Corporate Office:

F-196, Chhoriya Market, Raver, Dist: Jalgaon,
Maharashtra, India (PIN Code: 425508)

